

# Neural Lineage Markers Antibody Assembly kit

Neural Lineage Markers Antibody Assembly kit

Catalog # AP94038

## Product Information

<b>Application</b>	IHC-P, IF
<b>Host</b>	Rabbit
<b>Clonality</b>	Polyclonal


## Additional Information

<b>Dilution</b>	IHC-P=1:100-500,IF=1:100-500
<b>Format</b>	0.01M TBS(pH7.4) with 1% BSA, 0.09% (W/V) sodium azide and 50% Glyce
<b>Storage</b>	Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.

## Background


This product as supplied is intended for research use only, not for use in human, therapeutic or diagnostic applications.

## Images

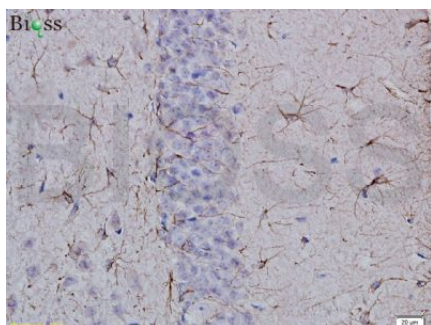


Primary: Anti- GFAP (bs-0199R) Tissue/cell: rat brain tissue; 4% Paraformaldehyde-fixed and paraffin-embedded; Antigen retrieval: citrate buffer ( 0.01M, pH 6.0 ), Boiling bathing for 15min; Block endogenous peroxidase by 3% Hydrogen peroxide for 30min; Blocking buffer (normal goat serum,C-0005) at 37°C for 20 min; Incubation: Anti-GFAP Polyclonal Antibody, Unconjugated(bs-0199R) 1:400, overnight at 4°C, followed by conjugation to the secondary antibody(SP-0023) and DAB(C-0010) staining

Primary: Anti-TUBB3 (Neuronal Marker) (bs-4512R) Tissue/cell: rat brain tissue; 4% Paraformaldehyde-fixed and paraffin-embedded; Antigen retrieval: citrate buffer ( 0.01M, pH 6.0 ), Boiling bathing for 15min; Block endogenous peroxidase by 3% Hydrogen peroxide for 30min; Blocking buffer (normal goat serum,C-0005) at 37°C for 20 min; Incubation: Anti-TUBB3/beta III Tubulin(Neuronal Marker) Polyclonal Antibody,

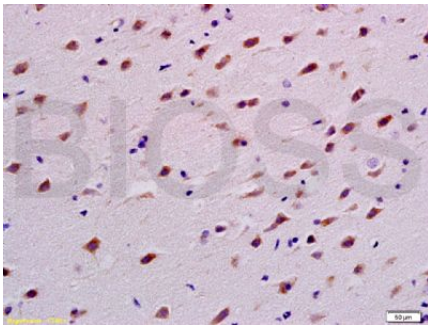


Unconjugated(bs-4512R) 1:200, overnight at 4°C, followed by conjugation to the secondary antibody(SP-0023) and DAB(C-0010) staining

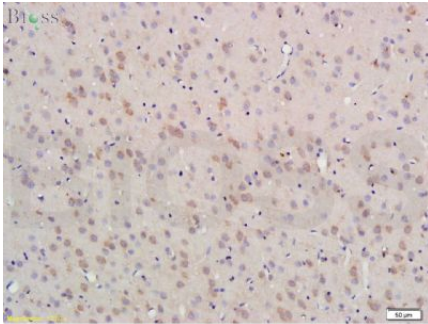


Primary: Anti-Nestin Polyclonal Antibody, Unconjugated(bs-0008R)


**bs-0008R Anti-Nestin**  
Formalin fixed and paraffin-embedded rat brain tissue labeled with Rabbit Anti-Nestin Polyclonal Antibody, Unconjugated(bs-0008R) at 1:300 followed by conjugation to the secondary antibody and DAB staining

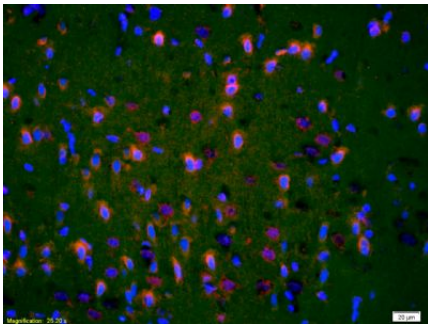

Primary: Anti-Nestin Polyclonal Antibody, Unconjugated(bs-0008R) Tissue/cell: rat brain tissue;4% Paraformaldehyde-fixed and paraffin-embedded; Antigen retrieval: citrate buffer ( 0.01M, pH 6.0 ), Boiling bathing for 15min; Blocking buffer (normal goat serum,C-0005) at 37°C for 20 min; Incubation: Anti-Nestin Polyclonal Antibody, Unconjugated(bs-0008R) 1:200, overnight at 4°C; The secondary antibody was Goat Anti-Rabbit IgG, Cy3 conjugated(bs-0295G-Cy3)used at 1:200 dilution for 40 minutes at 37°C. DAPI(5ug/ml,blue,C-0033) was used to stain the cell nuclei

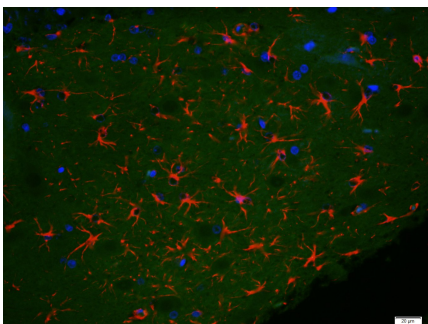


Primary: Anti- GFAP (bs-0199R) Paraformaldehyde-fixed, paraffin embedded (rat brain); Antigen retrieval by boiling in sodium citrate buffer (pH6) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (GFAP) Polyclonal Antibody, Unconjugated (bs-0199R) at 1:200 overnight at 4°C, followed by a conjugated secondary (bs-0295G-Cy3) at [1:500] for 90 minutes and DAPI staining of the nuclei.

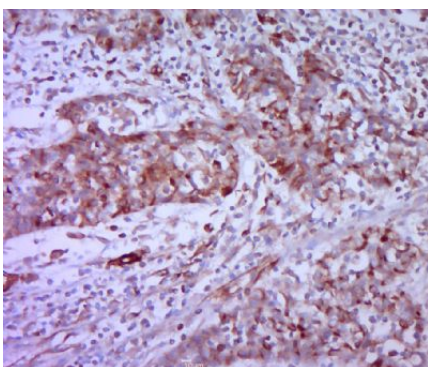

Primary: Anti- Vimentin (bs-0756R) Paraformaldehyde-fixed, paraffin embedded (human cervical cancer); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (Vimentin) Polyclonal Antibody, Unconjugated (bs-0756R) at 1:400 overnight at 4°C, followed by a conjugated secondary (sp-0023) for 20 minutes and DAB staining.

Primary:Anti-NeuN(bs-1613R) Paraformaldehyde-fixed, paraffin embedded (Mouse brain); Antigen retrieval by

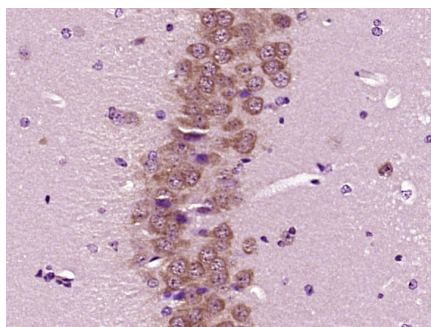

boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (NeuN) Polyclonal Antibody, Unconjugated (bs-1613R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.

Please note: All products are 'FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC OR THERAPEUTIC PROCEDURES'.