

# ox-LDL Rabbit pAb

ox-LDL Rabbit pAb Catalog # AP94170

#### **Product Information**

**Application** IHC-P, IHC-F, IF, ICC

Reactivity Mouse
Host Rabbit
Clonality Polyclonal
Calculated MW 31 KDa
Physical State Liquid

**Immunogen** Mouse Cu2+ oxidised low density lipoprotein

**Isotype** IgG

**Purity** affinity purified by Protein A

**Buffer** 0.01M TBS (pH7.4) with 1% BSA, 0.02% Proclin300 and 50% Glycerol.

**SUBCELLULAR LOCATION** Secreted.

**Important Note** This product as supplied is intended for research use only, not for use in

human, therapeutic or diagnostic applications.

**Background Descriptions** Low-density lipoprotein (LDL) is the carrier protein for cholesterol in the

blood. LDL binds to its receptor on the capillary walls and thereby mediates

the uptake and clearence of cholesterol from the circulation. In

atherosclerotic lesions oxidatively modified LDL is found and oxidized LDL is specifically recognized and ingested by macrophages via scavenger receptor A and CD36. Oxidized LDL may be a marker of atherosclerosis but the precise changes in oxidized LDL are not well described. Low-density lipoprotein

oxidised with Cu2SO4.

#### **Additional Information**

**Dilution** IHC-P=1:100-500,IHC-F=1:100-500,ICC/IF=1:100,IF=1:100-500

Format 0.01M TBS(pH7.4) with 1% BSA, 0.09% (W/V) sodium azide and 50% Glyce

**Storage** Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When


reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody

is stable for at least two weeks at 2-4 °C.

## **Background**

This product as supplied is intended for research use only, not for use in human, therapeutic or diagnostic applications.


### **Images**

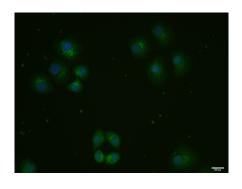

HepG2 cell; 4% Paraformaldehyde-fixed; Triton X-100 at room temperature for 20 min; Blocking buffer (normal goat serum, C-0005) at 37°C for 20 min; Antibody incubation with (ox-LDL) polyclonal Antibody, Unconjugated (AP94170) 1:100, 90 minutes at 37°C; followed by a conjugated Goat Anti-Rabbit IgG antibody at 37°C for 90 minutes, DAPI (blue, C02-04002) was used to stain the cell nuclei.



Paraformaldehyde-fixed, paraffin embedded Human Liver Cancer; Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15 min; Antibody incubation with ox-LDL Polyclonal Antibody, Unconjugated (AP94170) at 1:200 overnight at 4°C, followed by conjugation to the bs-0295G-HRP and DAB (C-0010) staining.

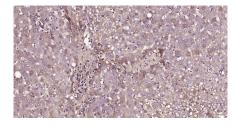

Paraformaldehyde-fixed, paraffin embedded Human Liver; Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15 min; Antibody incubation with ox-LDL Polyclonal Antibody, Unconjugated (AP94170) at 1:200 overnight at 4°C, followed by conjugation to the bs-0295G-HRP and DAB (C-0010) staining.

Please note: All products are 'FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC OR THERAPEUTIC PROCEDURES'.